



Edema
Management
Made Simple
with
FarrowWrap®

JOBST® FarrowWrap®



WHY JOBST® FarrowWrap®...

HOW IT WORKS

How does a JOBST® FarrowWrap® achieve compression?


- Utilizing special fabrics with smart elastic short-stretch bands comprised of multi-layered fabrics
- Optimal width of bands and amount of band overlap from bottom to top

It is important that every band of a JOBST® FarrowWrap® garment be applied at or near end-stretch. Having a consistent amount of tension throughout the garment helps to ensure graduated compression.

The Elastic Short-Stretch™ fabrics provide sustained resting compression while also augmenting the calf-muscle pump to deliver high working compression.


BENEFITS

- Easy use of the simple VELCRO® Brand hook & loop system saves therapists' and patients' time, and supports easy donning and doffing for maximum or optimal patient independence
- Since the product is adjustable, fluctuating edema levels can be more easily managed by the patients themselves, empowering them to control their own care
- The variety of styles and sizes available in Ready-to-Wear (RTW), Trim-to-Fit (TTF), Custom-Made (CM) ensures there is a FarrowWrap® to meet a broad range of clinical need

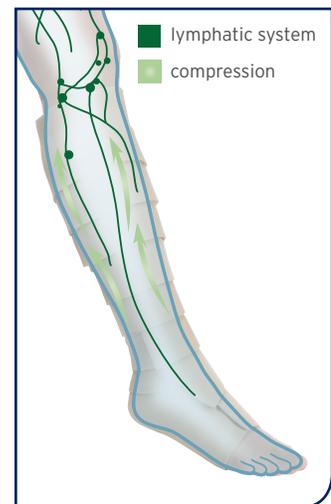


CHARACTERISTICS OF SHORT-STRETCH TECHNOLOGY

- Allows creation of products with different compression levels at end-stretch
- Provides augmentation of the body's natural muscle pump action
- Maximizes calf muscle function
- Helps minimize swelling throughout the day once garment is donned
- Provides higher compression when the patients stand or are active and less when they are not

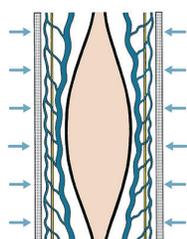

The pressure, exerted onto the body from the outside when the muscles at rest, is called "resting pressure"

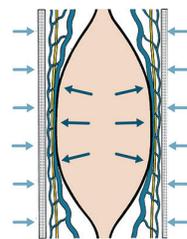

The pressure exerted onto the body from the inside as the muscle activates and expands, is called the "working pressure"

Short-stretch bandages deliver a low resting pressure and a high working pressure, producing a massaging effect on the calf muscle when the patient is active, which in turn produces a subsequent increase in lymphatic return. This makes short-stretch systems very effective in the management of lymphedema.

THE JOBST® FarrowWrap® PRODUCT RANGE

Adaptable and Affordable, Comfortable and Containing

Product Range


LITE

A double-laminated fabric suitable for mild to moderate compression needs, making it especially suitable for home care use / self-management patients. The fabric is both durable and breathable.

Available versions

- Upper and lower limb
- Ready-to-Wear (RTW)
- Trim-To-Fit (TTF)
- Custom-Made (CM)

20-30 mmHg* Mild to moderate chronic edema / lymphedema

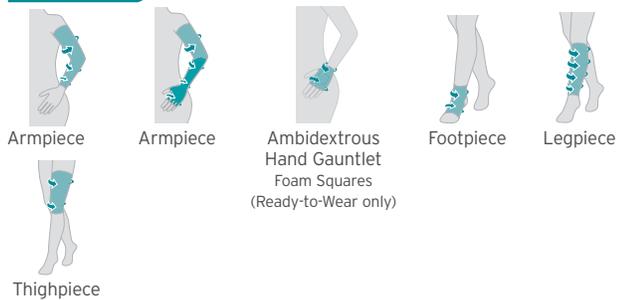

STRONG

Premium fabric with strong compression and soft inner layer suitable for patients with a good limb shape and moderate to severe edema.

Available versions

- Upper and lower limb
- Ready-to-Wear (RTW)
- Trim-To-Fit (TTF)
- Custom-Made (CM)

30-40 mmHg* Moderate to severe chronic edema / lymphedema with a good limb shape



CLASSIC

The original material that has been clinically proven over many years of use on thousands of patients worldwide.

Available versions

- Upper and lower limb
- Ready-to-Wear (RTW)
- Custom-Made (CM)

30-40 mmHg* Moderate to severe chronic edema / lymphedema



BASIC

A Value fabric with quality compression. Uses double-sided VELCRO® Brand product so it is possible to change band length and trim as limb reduces.

Available versions

- Lower limb
- Ready-to-Wear (RTW)

30-40 mmHg* Moderate to severe chronic edema / lymphedema



4000

The inner sleeve and 4-band design make for easy application. Anatomically contoured to prevent gapping.

Available versions

- Lower limb
- Ready-to-Wear (RTW)

30-40 mmHg* Moderate to severe chronic edema / lymphedema



PRODUCT ASSORTMENT

	Footpiece	Legpiece	Foot Compression Liner	Thighpiece	Armpiece	Ambidextrous Hand Gauntlet Foam Squares
Mild to Moderate 20-30 mmHg						
JOBST® FarrowWrap® LITE	RTW, CM, TTF	RTW, CM, TTF		RTW, CM	RTW, CM, TTF	RTW
JOBST® FarrowHybrid® ADI			RTW			
Moderate to Strong 30-40 mmHg						
JOBST® FarrowWrap® STRONG	RTW, CM, TTF	RTW, CM, TTF		RTW, CM	CM	
JOBST® FarrowWrap® CLASSIC	RTW, CM	RTW, CM		RTW, CM	CM	
JOBST® FarrowWrap® 4000		RTW				
JOBST® FarrowHybrid® ADII			RTW			
JOBST® FarrowWrap® BASIC	RTW	RTW				

RTW: Ready-To-Wear • CM: Custom-Made • TTF: Trim-To-Fit

Accessories


Hand Gauntlet Foam Squares to be used with JOBST® FarrowWrap® Lite Hand Gauntlet. The foam inserts help reduce edema that can pool on top of the hand.



JOBST® FarrowWrap® GarmentGrip™ can be cut to size to fit the shape and contour of the body and VELCRO® brand hook-and-loop fasteners can be added to secure in place. For triangular-shaped legs, it is recommended that it be worn over the legpiece to secure the wrap system in place.

Pad & Shape

JOBST® FarrowFoam comes in a variety of sizes and rolls or sheets. It is great for shaping any irregular limb shapes or providing padding to sensitive areas.

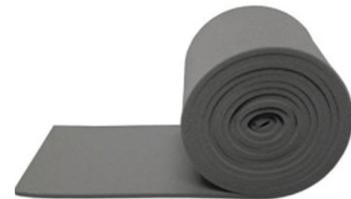

Ready-To-Wear (RTW)

Garments are pre-made pieces to fit the majority of patient sizes.

Trim-To-Fit (TTF)

Garments offer one garment for both active/reduction and maintenance therapy. A Trim-To-Fit garment can be trimmed to adapt as extremity reduces. It can fit most patients in a single size.

Custom-Made (CM)

Garments are manufactured to the specific limb measurements of each patient to provide a better fit for those limb shapes that do not fit within the RTW size range.

UPPER EXTREMITY

JOBST® FarrowWrap® ARMPIECE - SUITABLE FOR 24-HOUR WEAR

LITE RTW Armpiece*



LITE RTW armpiece is an anatomically contoured, effective solution to control upper extremity edema. It is easy to apply and remove with one hand and designed for mild to moderate edema control.

- 20-30 mmHg

LITE TTF Armpiece*



The LITE TTF armpiece can be customized at the point of purchase and will fit most limb sizes. The LITE TTF armpiece can be retrimmed as the patient's edema reduces.

- 20-30 mmHg

LITE Custom Armpiece*



The LITE custom armpiece consists of overlapping bands of short-stretch fabric. It generally requires two hands to apply. The short-stretch bands allow garment to be applied at end stretch.

- 20-30 mmHg

CLASSIC Custom Armpiece*



CLASSIC custom armpiece is for more severe edema requiring high levels of compression. It is only available as a custom garment and generally requires two hands to apply. Patients with very large arms or moderate to severe swelling would benefit from a CLASSIC garment.

- 30-40 mmHg

STRONG Custom Armpiece*



The STRONG custom armpiece is used in the same situations as CLASSIC but has the added benefit of ease of maintenance.

- 30-40 mmHg

* Colours: Tan, Black; Sizes: S, M, L; Available in Left or Right, and Short, Regular & Long Lengths

JOBST® FarrowWrap® HANDPIECE - SUITABLE FOR 24-HOUR WEAR



JOBST® FarrowWrap® LITE Hand Gauntlet 20-30mmHg*

JOBST® FarrowWrap® LITE Hand Gauntlet is ambidextrous and comes with 4mm and 8mm foam squares for added compression to the dorsal hand. Suitable for mild to moderate edema.

Benefits that give your patient greater freedom to enjoy life:

- Foam inserts help reduce edema that can pool on the top of the hand
- Ambidextrous, i.e. can be used on either hand
- Durable and comfortable

LINER



Flesh Tone Liner**

Liner is only required for trim-to-fit custom armpiece. There is no compression in liner.

The liner provides a smooth base and skin protection prior to applying the outer wrap piece. In addition, daily washing of the outer wrap piece is not recommended due to the nature of the material and the VELCRO® brand hook-and-loop fasteners. However, a fresh liner can be worn daily as the liner can be washed as frequently as needed.

- Contains silver

* Colours: Tan; Sizes: XS, S, M, L

** Colours: Tan; Sizes: S, L

LOWER EXTREMITY

JOBST® FarrowWrap® THIGHPIECE AND LEGPIECE - SUITABLE FOR 24-HOUR WEAR

LITE Legpiece**



CLASSIC Legpiece**



STRONG Legpiece**



BASIC Legpiece**



LITE Thighpiece



CLASSIC Thighpiece



STRONG Thighpiece



JOBST® FarrowWrap® LITE* is designed for mild to moderate edema control. Available in Ready-To-Wear (RTW), Trim-To-Fit (TTF), and Custom-Made (CM). LITE garments can be machine washed and dried.†

- 20-30 mmHg

JOBST® FarrowWrap® CLASSIC* is designed for moderate to severe edema control. It is available in Ready-To-Wear (RTW) and Custom-Made (CM) garments.

Hand wash only.†

- 30-40 mmHg

JOBST® FarrowWrap® STRONG* is designed for moderate to severe edema control and is available in Ready-To-Wear (RTW) and Custom-Made (CM) sizes. STRONG garments can be machine washed and dried.†

- 30-40 mmHg

JOBST® FarrowWrap® BASIC* is an exceptional value and features trimmable bands and adjustable VELCRO® Brand products so the garment can be trimmed as the limb reduces. BASIC garments can be machine washed and dried.†

- 30-40 mmHg

Available in Legpiece only.

* Colours: Tan; Sizes: XS, S, M, L, XL; Available in Short, Regular and Tall Lengths

** Colours: Tan; Sizes: XS, S, M, L, XL; Available in Regular and Tall Lengths

† Follow care instructions in the product packaging.

JOBST® FarrowWrap® LEGPIECE - SUITABLE FOR 24-HOUR WEAR



JOBST® FarrowWrap® 4000

JOBST® FarrowWrap® 4000 legpiece is designed for moderate to severe edema and is easy to apply with just four bands and inner sleeve. It is available in RTW sizing. 4000 garments can be machine washed and dried.



Flesh Tone Liner**

- Fits up to 70cm calf circumference
- Contains silver

The liner provides a smooth base and skin protection prior to applying the outer piece, and can be washed daily.

Customizing garments made easy with Trim-to-Fit

TTF garments can be customized on site to fit a variety of patients. You can easily retrim the garment if the patient's edema reduces. There are two TTF sizes which cover 6 standard sizes. TTF is available in STRONG and LITE fabrics.

On triangular shaped legs we recommend **JOBST® FarrowWrap® GarmentGrip™** be worn over our legpiece to help hold the garment in place.

* Colours: Tan, Black; Sizes: XS, S, M, L
** Colours: Tan; Sizes: S, L

LOWER EXTREMITY

JOBST® FarrowWrap® FOOTPIECE - SUITABLE FOR 24-HOUR WEAR

LITE Footpiece*



LITE footpiece provides compression for mild to moderate edema.

- 20-30 mmHg

CLASSIC Footpiece*



CLASSIC footpiece is clinically proven in managing edema. Durable construction.

- 30-40 mmHg

STRONG Footpiece*



STRONG footpiece provides compression for moderate to severe edema.

- 30-40 mmHg

BASIC Footpiece*



BASIC footpieces provide compression for moderate to severe edema.

- 30-40 mmHg

Customizing garments

TTF footpieces can be instantly customised using our proprietary TTF innovation. Must be trimmed before donning.

* Colours: Tan; Sizes: XS, S, M, L, XL; Available in Regular and Long Lengths

THE THERAPEUTIC CONTINUUM OF CARE FOR EVERY PATIENT JOURNEY

It is critical to take proactive measures using an integrated approach combining compression therapy with treatments that support optimal wound healing. Essity supports the Continuum of Care approach that enhances outcomes by combining therapies and choosing patient-centered devices to optimize results.



REMOVE EXUDATE & BACTERIA



Cutimed® Sorbion® Sorbact®
featuring Sorbact® Technology



Cutimed® Sorbion Sachet®
Superabsorbent Dressing



COMPRESS & REDUCE EDEMA



JOBST® Compri 2
Two-Layer Bandaging

JOBST® Elvarex®
Custom Compression
Garments



JOBST® Opaque
Ready-to-Wear Hosiery



JOBST® FarrowWrap®
Adjustable Short Stretch
Compression

To learn more about
JOBST® products and services, call
1-877-978-5526
or visit
www.jobstcanada.com | www.essity.com

* The average compression for a medium sized ankle.
VELCRO® is a Registered Trademark of Velcro BVBA.
Please refer to the product label and / or package insert for full instructions on the safe use of these products.



JOBST®,
an Essity brand

132666 K9 - 0.75M

Essity
1275 North Service Road West, Suite 800,
Oakville, ON L6M 3G4

www.jobstcanada.com
www.essity.com
Tel 1-877-978-5526
Fax 1-877-978-9703